
Product model 8100A

product structure

1. massage head
 2. indirect vibration working mode indicator
 3. heating working indicator
 4. vibratory working mode indicator
 5. switch on
 6. USB charging cable
 7. USB charging port
- Please use 5v/1A adapter to recharge

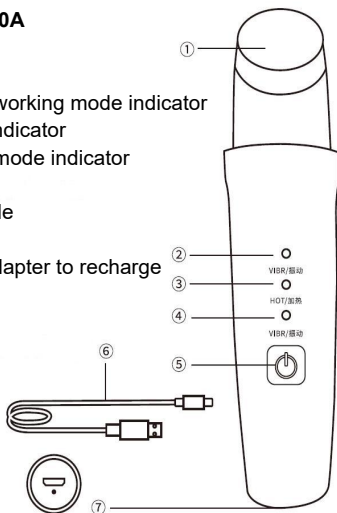Product introduction

1. Thermal beauty eye and instrument with constant temperature 42degree and high-frequency vibration principle to introduced the skin care products, nutrition into the skin bottom layer,play an efficient absorption role, so as to achieve eye,lip,facial beauty effect.
- 2.Small appearance, easy to carry,US Bcharging,safe and environmentally friendly
- 3.Product massage head use special big head design which can fie every angle of skin,with high frequency vibration,constant temperature 42c, indirect vibration,high-frequency vibration+constant temperature and indirect vibration+constant temperature all five mode can be used for lip lines, eye dark circles.

Eye care+lip care:

- 1.apply eye skin care products(such as eye cream,eye mask,etc) evenly around the eyes
- 2.apply lip care products(such as lipstick,eye serum,etc)evenly to the lips.
- 3.long press the key for 3 seconds, the VIBR indicator light turns on blue light vibration mode to start,then short press on the key once,the functional conversion indicator lights up to blue,HOT constant temperature mode start,then short press the key once,VIBR indicator light blue indirect vibration mode start,again short press on the key one time function conversion VIBRA+HOT indicator light bright blue light,high frequency vibration+constant temperature mode start:then short press the key once function conversion VBR+HOT indicator light blue ,indirect vibration mode+constant heat mode start.
- 4.after choose the mode,massage head from the eyes to the temple gently point and press to promote absorption
- 5.the product stops automatically after 5 minutes in each mode
- 6.the nursing mode can also be changed according to individual needs.

NOTE ITEMS

- 1.before using the product, please read the instructions carefully and keep them for future reference
- 2.the product is charged with USB cable, it takes 1 hour to charge when the battery in low.
- 3.when charging, the 3 lights turn blue and flash in sequence.When full ,3 indicator lights turn blue, no longer flashing.
- 4.If you feel uncomfortable during using.Or you skin appears abnormal, please stop using the well and consult your doctor immediately. Products without water function, do not place in water.
- 5.please put the product in a dry and ventilated place
- 6.avoid direct sunlight no corrosive chemicals around.
- 7.when using,please take proper amount of skin care products,the skin can not absorb all leading to waste at the same time.There is the possibility of residual infiltrating into the product, damage the product life.Please clean it immediately after use wipe the cloth clean.
- 8.do not let children touch the product in order to avoid danger
- 9.the surface of the appliance is hot, and people who are not sensitive to heat must pay attention to it when using it.
- 10.before discarding the appliance , the battery must be removed from the appliance: the appliance must be powered off when the battery is removed .batteries should be disposed of safely.

- 1.Press 3 seconds to turn on(Note:need to unplug the power cord)
- 2.the machine about 5 minutes chip will automatically shutdown, need to continue to use can restart
- 3.charge time for about 2 hours,indicator lights three long do not flash indicating full

3Pagen Versand und Handelsges. mbH
St. Jöris Strasse 16-28
52477 Alsdorf
kontakt@3pagen.de

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