Rechargeable Heated Hand Massager

Item No: P655 / 624850

Instructions for Use

v001: 11/03/25

Please read and retain these instructions for future reference

Java

wellness

ENGLISH INSTRUCTIONS

Power Details:

Voltage details: 100-240V~

Frequency details: 50Hz-60Hz

Wattage details: 6W

Accessories:

1x hand massager

1x adaptor

1x cable

1x manual

Please read this instruction manual thoroughly before starting and ensure you are familiar with the operation of your new RECHARGABLE HEATED HAND MASSAGER.

User Instructions here

1. Control panel area
(Time setting/Pressure/Heat/Decompression/Mode)
2. Digital tube display screen
3. DC jack (charging port)


Charge operation:

1. This item works with a rechargeable battery and to charge, please insert the cylinder plug of DC charging line into the DC power supply socket of the product. The other end is connected to 5V/1A-4A charger or 5V/1A adapter.
2. While charging the indicator light will appear red and once fully charged the light will appear blue.
3. The Battery can be used for about 2-3 hours (without heat) after full charge.
4. This item can be used while charging.

Use operation:

1. Spread your fingers and palms down into the inner gallbladder of the massager.
2. Turn on the device
3. Select one of the operation keys and enjoy the massage.

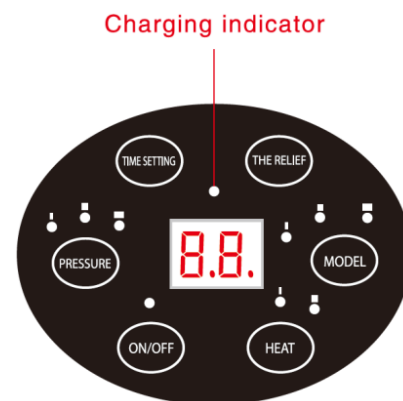Pressure Key: Controlling massage intensity, this machine has three kinds of massage intensity, which is displayed by indicator lamp. (I is low gear, II is medium gear and III is high gear).

Heating key: The heating can be controlled independently, when the massage function is turned on. There are two heat functions (I is low heat, II is medium heat)

Decompression key: Controls the decompression function of the whole machine. One short press of this button and the machine will decompress the air bags and free your hand.

Mode key: Control mode, there are three kinds of massage modes in this machine (1. health care mode, 2. vitality mode and 3. relaxation mode).

4. **Time Setting:** Time setting button lets you select a 5-, 10- or 15-minute massage. Press button once for 5 minutes, twice for 10 minutes and three times for 15 minutes. The device will automatically stop once the time is up.

Attention:

1. Do not use hand massager while driving.
2. Do not use massager in high temperature and humidity.
3. Do not use If you have a hand/ finger injury or have recently undergone surgery.
4. Make sure your palms are not Jewelled before use.
5. We recommend not to apply nail polish an hour before use.
6. It is recommended to wear disposable gloves to use this product.
7. When the storage period exceeds 6 months, it should be checked and power-on.
8. Do not use if you feel unwell.
9. Hand massager is designed for household electronic health care products. It is not suitable for clinics, treatment and other commercial purposes.
10. If cleaning is required, sterile cotton swab soaked with neutral detergent should be used. No Cleaning Products with Water.
11. This massager has heating function. Please pay attention to users who are insensitive to temperature when using it.

WARNINGS AND SAFETY INSTRUCTIONS:

WHEN USING ELECTRICAL APPLIANCES, BASIC SAFETY PRECAUTIONS SHOULD ALWAYS BE FOLLOWED:

- READ ALL INSTRUCTIONS BEFORE USE.
- FOR DOMESTIC INDOOR USE ONLY.
- DO NOT USE OUTDOORS OR FOR COMMERCIAL USE.
- ONLY USE WITH THE MAINS ADAPTOR/CHARGER SUPPLIED.
- **DO NOT** USE ANY OTHER MAINS ADAPTOR/CHARGER AS THIS MAY CAUSE DAMAGE AND INVALIDATES ANY WARRANTY YOU MAY HAVE.
- BEFORE USE, CHECK THOROUGHLY FOR ANY DEFECTS AND DO NOT USE IF DEFECTS ARE FOUND. TAKE CARE NOT TO DROP THE APPLIANCE AS HEAVY IMPACTS MAY CAUSE INTERNAL DAMAGE.
- BEFORE CONNECTING THE APPLIANCE, CHECK IF THE VOLTAGE INDICATED ON THE BOTTOM OF THE APPLIANCE IS IN ACCORDANCE WITH THE MAINS VOLTAGE IN YOUR HOME.
- CHECK THE APPLIANCE AND POWER CABLE CAREFULLY FOR DAMAGE BEFORE EACH USE.
- THIS APPLIANCE CAN BE USED BY CHILDREN AGED FROM 8 YEARS AND ABOVE, AND PERSONS WITH REDUCED PHYSICAL, SENSORY OR MENTAL CAPABILITIES OR LACK OF EXPERIENCE AND KNOWLEDGE IF THEY HAVE BEEN

GIVEN SUPERVISION OR INSTRUCTION CONCERNING USE OF THE APPLIANCE IN A SAFE WAY AND UNDERSTAND THE HAZARDS INVOLVED. CHILDREN SHALL NOT PLAY WITH THE APPLIANCE. CLEANING AND USER MAINTENANCE SHALL NOT BE MADE BY CHILDREN WITHOUT SUPERVISION.

- CHILDREN OF LESS THAN 3 YEARS SHOULD BE KEPT AWAY UNLESS CONTINUOUSLY SUPERVISED.
- **SOME PARTS OF THIS PRODUCT CAN BECOME VERY HOT AND CAUSE BURNS. PARTICULAR ATTENTION HAS TO BE GIVEN WHERE CHILDREN AND VULNERABLE PEOPLE ARE PRESENT.**
- THE USE OF ANY ACCESSORY OR ATTACHMENTS NOT RECOMMENDED BY COOPERS OF STORTFORD MAY CAUSE INJURIES AND INVALIDATE ANY WARRANTY YOU MAY HAVE.
- DO NOT LET THE POWER CABLE HANG OVER THE EDGE OF A TABLE OR COUNTER, OR TOUCH HOT SURFACES.
- DO NOT PLACE ON OR NEAR HEAT SOURCES.
- USE ON A LEVEL, STABLE SURFACE.
- DO NOT COVER THE ITEM IN ANY WAY WHEN IN USE AS THIS MAY CAUSE OVERHEATING.
- DO NOT USE THIS APPLIANCE FOR ANY OTHER USE THAN THE INTENDED USE IT IS DESIGNED FOR.
- ENSURE THERE IS SUFFICIENT DISTANCE BETWEEN THIS APPLIANCE AND OTHER ITEMS ON THE WORK SURFACE OR WORK STATION, OR SUFFICIENT DISTANCE BETWEEN THE APPLIANCE AND THE CONFINES OF THE AREA IN WHICH IT IS SITUATED. WE RECOMMEND A MINIMUM CLEARANCE AROUND ALL FACES OF AT LEAST 10CM.
- KEEP ITEM AT LEAST 1M AWAY FROM COMBUSTIBLE MATERIAL SUCH AS FURNITURE, CUSHIONS, BEDDING, PAPER, CLOTHES, CURTAINS ETC.
- TO PROTECT AGAINST ELECTRIC SHOCK, DO NOT IMMERSE POWER CABLE, MAINS ADAPTOR/ CHARGER OR APPLIANCE, IN WATER OR ANY OTHER LIQUID.
- UNPLUG FROM OUTLET WHEN NOT IN USE AND BEFORE CLEANING.
- DO NOT OPERATE ANY APPLIANCE WITH A DAMAGED POWER CABLE OR MAINS ADAPTOR/ CHARGER OR AFTER THE APPLIANCE MALFUNCTIONS OR HAS BEEN DAMAGED IN ANY MANNER. WE RECOMMEND THAT THE POWER CABLE IS NOT KINKED OR STORED IN A TWISTED MANNER, AND THAT YOU REGULARLY CHECK ITS CONDITION.
- DO NOT LEAVE UNATTENDED.
- ALLOW TO COOL BEFORE CARRYING OR PLACING IN STORAGE.
- ANY REPAIRS MUST BE CONDUCTED BY A QUALIFIED ELECTRICIAN OR ELECTRICAL SERVICE AGENT.

PLUG WIRING (UK & IRELAND)

This appliance is supplied with a mains adaptor/charger with BS1363 plug pins.

DO NOT use if the socket outlets in your home or office are not suitable for this product's mains adaptor/ charger plug pins. Attempts to insert the mains adaptor/charger into the wrong socket is likely to cause electric hazard or damage.

DO NOT use with any other mains transformer/charger than the one that is supplied with this product.

CARE AND CLEANING INSTRUCTIONS

- Always switch off and unplug before cleaning.
- Do NOT immerse in water. Risk of electric shock.
- Do NOT use abrasive cleaners or sharp utensils to clean this appliance.
- Wipe surfaces with a soft dry cloth.
- Allow to dry thoroughly before further use.

DISPOSAL

- Coopers of Stortford use recyclable or recycled packaging where possible.
- Please dispose of all packaging, paper, cartons in accordance with your local authority recycling regulations.
- At the end of the product's lifespan please check with your local council authorised household waste recycling centre for disposal.

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The Enterprise Department Ltd, Trading as Coopers of Stortford.
The Enterprise Department Ltd, Trading as -